

Customised Monitor Kit for EP50 Camera

See Your Sample Through a Full HD Display



SZ Series

SZX Series

CKX Series

CX Series



Tripod Adapter

Article Number: E9701278

12 inch Full HD Screen Monitor Kit

Article Number: E0437042

Discover the Benefits

- **See in HD** – view high resolution images and videos in a full HD display
- **Collaborate in comfort** – say goodbye to uncomfortable positions associated with ocular lens viewing with on-screen imaging
- **Learn together** – improve student learning experience by enabling side-by-side viewing
- **Safety first** – wear full safety protection without compromising visualization
- **Set-up made simple** – set-up your high definition monitor in 5 simple steps

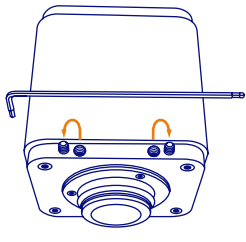
Customised Monitor Kit for EP50 Camera

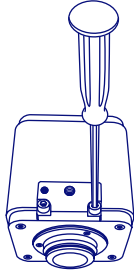
Key Customer Applications

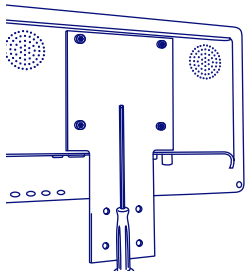
Spanning multiple applications – including material inspection, cytology, hematology, and cell culture – the EP50 camera makes digital imaging of microscopy samples a fast and engaging experience. Now, you can couple your EP50 camera with a high quality monitor to take your application and your comfort to the next level.

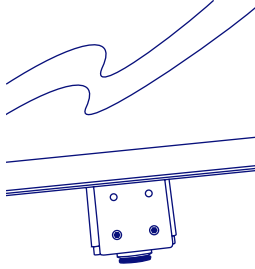
- Cytology, Hematology
- Cell Culture
- Education / Training in Course rooms in Medical Schools and Universities
- Material Inspection in Industry
- Sample Processing Under the Stereo Microscope

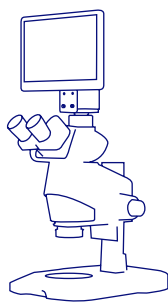
5 Simple Steps to Set Up the Monitor

- 

1. Remove the headless screws from the EP50
- 

2. Attach the tripod adapter on to EP50
- 

3. Attach adapter on to the monitor
- 

4. Attach monitor adapter to the tripod adapter on the EP50
- 

5. Install the setup on to a microscope

Single Kit That Contains All Required Items for Multiple Configurations



Standard Configuration

Directly attach your monitor to the top of the microscope using the included short HDMI cable for:

- A small and compact footprint
- Ergonomic observation



Benchtop Stand Configuration

Attach your monitor using the benchtop stand with the included long HDMI cable to:

- Extend the display
- Give a better view of the stage
- Allow multiple observers



 www.olympus-europa.com

For research only. This product is not intended to be used for IVD.
Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.

OLYMPUS

OLYMPUS EUROPA SE & CO. KG

Postbox 10 49 08, 20034 Hamburg, Germany
Amsinckstrasse 63, 20097 Hamburg, Germany
Phone: +49 40 23773-0, Fax: +49 40 233765
E-Mail: scientificolutions@olympus-europa.com